


Part E

1. Dial pen to **your** insulin dose.

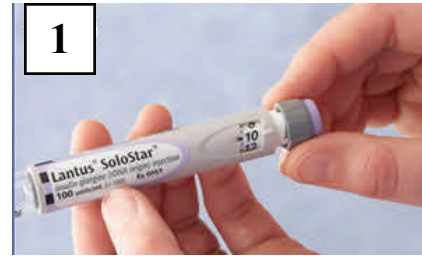

Part F

1. Inject needle into the site you have chosen.
2. Press button on bottom of pen until you hear it click and the dose in the dose window moves back down to zero.
3. Leave needle in your skin and keep pressing the button for at least 5-10 seconds .
4. Remove needle from your skin.



Part G

1. Put large cap back on needle and unscrew counter-clockwise to remove pen needle.
2. Discard needle in hard container (for example: sharps container, empty detergent or shampoo bottle). (When container is full, tape lid and discard in trash.)
3. Put cap back on pen.



Part H

1. Keep the pen that you are using at room temperature.
It is good for:
 - 56 days if Tresiba
 - 42 days if Levemir or Toujeo
 - 28 days if Lantus, Basaglar, Humalog, Novolog or Apidra
 - 2 weeks if NPH
 - 10-14 days if a mixed insulin, Novolog 70/30, or Humalog 75/25

Keep pens that you have not started using in the refrigerator.